


Question or Item	Answer
Did you verify that you are using the most recent version of this form?	<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No
Date	2024-10-08
Individual creating the CSB	Brad Hunt
Customer:	All
Product(s) affected:	Welch Allyn Connex® Spot Monitor (CSM) Welch Allyn Spot Vital Signs® Monitor (Spot4400)
Hardware Version(s) affected:	ALL
Software Version(s) affected:	ALL
Subject:	FRIWO Power Supply: Adapter Plate Removal

Release Details

Starting in October 2024, the CSM and Spot 4400 Devices will no longer be shipped with the 421008 Power Supply Adapter Plate and the 4400-MBS and 7000-MS3 mobile stands will no longer be shipped with the 411626 Power Supply Adapter Plate (Pictured Below):

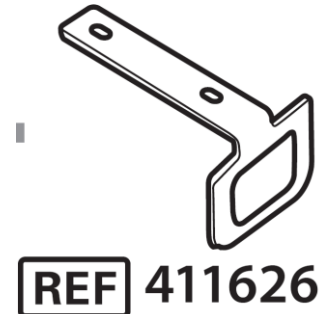

The 421008 adapter plate was introduced to the devices in March 2022 to allow backwards compatibility of the FRIWO device power supply with legacy mobile stand designs.



The 411626 adapter plate was introduced to the mobile stands in March 2022 to allow backwards compatibility of the mobile stands with legacy FRIWO device power supply designs.



The 421008 and 411626 adapter plates will now be available for individual sale as service kits. For assembly instructions of the 421008 and 411626 adapter plates, please see the CSM and/or Spot 4400 Service Manual.

END OF BULLETIN

Version	Sec, Pg, Para Changed	Change Made	Date Version Created	Version Created By (initials)
A	N/A	Initial Release	2024-10-08	BH